

# Large Capacity Refrigerated Centrifuge

## Union 32R




**Union 32R**  
**10,000rpm**  
**13.975 xg**  
**4x 750ml**  
**-20°C to +40°C**

This has the same functions as the Union 5KR except for a smaller capacity (104 samples per each operation). However, the Union 32R is so compact that you can use it on your laboratory table.

### Features

- one chip microprocessor control
- Maintenance-free, quiet, Brushless AC induction motor
- 10 program storage and recall
- Display and control RCF
- Gentle start-up and braking : ten acceleration/deceleration profiles
- RCF/RPM calculation
- Rotor recognition
- Temperature limit input to ensure proper temperature control
- Display abnormality : imbalance, over-speed, over-temperature, door-open

### S750-4B

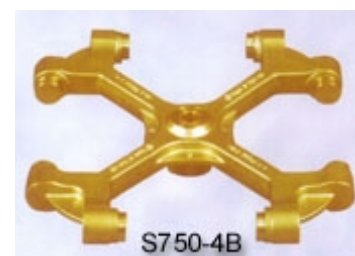

### Specification

<b>Max. RPM</b>	10,000 rpm
<b>Max. RCF</b>	13,975 xg
<b>Max. Capacity</b>	4 x 750ml
<b>Speed preset and display</b>	1 rpm
<b>Temperature Range</b>	-20°C to +40°C
<b>Temperature limit set range</b>	0°C to +/- 30°C
<b>Timer</b>	9hurs 59min 59sec and "Hold" for free run
<b>Deceleration time</b>	10 setting
<b>Programmability</b>	10 memory
<b>Rotor Recognition</b>	10 rotor
<b>Digital display</b>	RCF,RPM, Time. Program, Rotor Number, Temperature Temp. Limit, Brake
<b>Power</b>	230Vac, 50/60Hz, Single phase, 1.5 Kw
<b>Dimension</b>	640W x 630D x 450 H/mm

*labkorea*