


PRECISION INCUBATOR



PI-80H


- Well designed air mixing system provides accurate temperature control.
- Digital temperature control.
- Rubber magnet seal door.
- Powerful and quiet DC motor.
- Compact Rocker(CR95) can be operated inside the incubator with the electrical socket provided.
- Movable sensor :Thermo sensing probe can be closely placed to samples on Rocker platform.

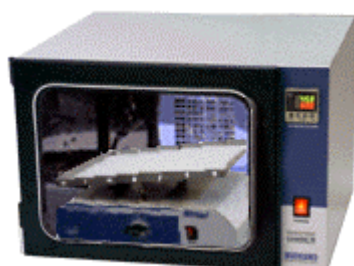

SPECIFICATION

SPECIFICATION	
MODEL	PI-80H
Temperature range	Ambient +5℃ ~ 80℃
Temperature accuracy	± 0.1℃
Dimension(mm)	out side: W470 × D418 × H310 (W18.5" × D16.5" × H12.2") in side: W 336 × D250 × H250 (W13.2" × D9.8" × H9.8")
Weight	19kgs(42 Lbs)
Power	AC 110V (AC 220V)

ACCESSORIES

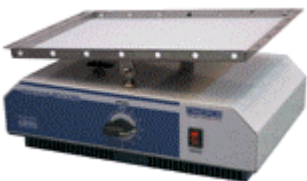Accessory	model	Description
	CR95	Able to use the CR95, compact ROCKER, inside PI-80H

labkorea